

VetLine Brucella

ELISA

Enzyme immunoassay for the qualitative determination of antibodies against Brucella in veterinary serum, pooled serum samples and milk samples.

Enzymimmunoassay zur qualitativen immunenzymatischen Bestimmung von Antikörpern gegen Brucella in veterinärem Serum, gepoolten Serumproben und Milchproben.

Enzimoimmunoensayo para la determinación cualitativa de anticuerpos contra Brucella en suero veterinario, muestras del pool de suero y muestras de leche.

English	2
Deutsch	7
Español	12
Bibliography / Literatur / Bibliographie / Bibliografía / Bibliografia	18
Abbreviations / Abkürzungen / Abréviations / Abbreviazioni / Abreviaciones / Abreviaturas	18
Symbols Key / Symbolschlüssel / Explication des Symboles / Legenda / Símbolos / Tabela de símbolos	19
Summary of Test Procedure / Kurzanleitung Testdurchführung / Résumé de la procedure de test / Schema della procedura / Resumen de la técnica / Resumo do Procedimento de Teste	20

Product Number: BRUVT0050 (96 Determinations)

To the online shop:

ENGLISH

1. INTRODUCTION

Brucella species are very small, gram-negative immobile bacteria (0.4 - 0.8 microns in diameter and 0.4 - 3.0 microns in length). They were named after the English military doctor Bruce Davis, who isolated the pathogen 1887 on Malta from the spleen of a deceased soldier with undulating fever.

Especially following four species are of importance:

<i>Brucella abortus</i> :	causative agent of bovine brucellosis
<i>Brucella melitensis</i> :	causative agent of ovine and caprine brucellosis
<i>Brucella suis</i> :	causative agent of porcine brucellosis
<i>Brucella canis</i> :	causative agent of canine brucellosis

The symptoms of the disease are similar in cattle, sheep, goats and pigs. Sheep seem to be less susceptible so less frequently abortion is observed in this species. After infection, brucellosis begins with a clinically normal phase which is often unnoticed. Later there are often joint inflammation, and sometimes mastitis. The main symptoms of brucellosis are abortions, premature births and birth of dead or weak offspring. Infected bulls can develop orchitis and epididymitis. The pathogen is transmitted mainly by ingestion and mating. Infection by milk, urine and faeces is also possible. In livestock brucellosis can become epidemic.

Brucellosis is a transmitted disease to humans (zoonosis). Individuals in risk are especially shepherds, farmers, animal keepers, veterinarians and laboratory staff.

Infections may be diagnosed:

- By genus specific PCR
- Serology: Detection of antibodies by ELISA

2. INTENDED USE

The NovaTec VetLine *Brucella* ELISA is intended for the qualitative determination of antibodies against *Brucella* in veterinary serum, pooled serum samples and milk samples.

3. PRINCIPLE OF THE ASSAY

The qualitative immunoenzymatic determination of specific antibodies is based on the ELISA (Enzyme-linked Immunosorbent Assay) technique.

Microtiterplates are coated with specific antigens to bind corresponding antibodies of the sample. After washing the wells to remove all unbound sample material a horseradish peroxidase (HRP) labelled conjugate is added. This conjugate binds to the captured antibodies. In a second washing step unbound conjugate is removed. The immune complex formed by the bound conjugate is visualized by adding Tetramethylbenzidine (TMB) substrate which gives a blue reaction product.

The intensity of this product is proportional to the amount of specific antibodies in the sample. Sulphuric acid is added to stop the reaction. This produces a yellow endpoint colour. Absorbance at 450/620 nm is read using an ELISA Microtiterplate reader.

4. MATERIALS

4.1. Reagents supplied

- **Microtiterplate:** 12 break-apart 8-well snap-off strips coated with *Brucella* antigens; in resealable aluminium foil.
- **Sample Dilution Buffer:** 1 bottle containing 100 mL of phosphate buffer (10 mM) for sample dilution; pH 7.2 ± 0.2; coloured yellow; ready to use; white cap; ≤ 0.0015 % (v/v) CMIT/MIT (3:1).
- **Stop Solution:** 1 bottle containing 15 mL sulphuric acid, 0.2 mol/L; ready to use; red cap.
- **Washing Buffer (20x conc.):** 1 bottle containing 50 mL of a 20-fold concentrated phosphate buffer (0.2 M), pH 7.2 ± 0.2, for washing the wells; white cap.
- **Conjugate:** 1 bottle containing 20 mL of peroxidase labelled Protein A/G; coloured yellow; ready to use; white cap; ≤ 0.02 % (v/v) MIT.
- **TMB Substrate Solution:** 1 bottle containing 15 mL 3,3',5,5'-tetramethylbenzidine (TMB), < 0.1 %; ready to use; yellow cap.
- **Positive Control:** 1 vial containing 2 mL; coloured yellow; ready to use; red cap; ≤ 0.02 % (v/v) MIT.
- **Cut-off Control:** 1 vial containing 3 mL; coloured yellow; ready to use; green cap; ≤ 0.02 % (v/v) MIT.
- **Negative Control:** 1 vial containing 2 mL; coloured yellow; ready to use; blue cap; ≤ 0.0015 % (v/v) CMIT/MIT (3:1).

For hazard and precautionary statements see 12.1

For potential hazardous substances please check the safety data sheet.

4.2. Materials supplied

- 1 Cover foil
- 1 Instruction for use (IFU)
- 1 Plate layout

To the online shop:



4.3. Materials and Equipment needed

- ELISA Microtiterplate reader, equipped for the measurement of absorbance at 450/620 nm
- Incubator 37 °C
- Manual or automatic equipment for rinsing Microtiterplate wells
- Pipettes to deliver volumes between 10 and 1000 µL
- Vortex tube mixer
- Distilled water
- Disposable tubes

5. STABILITY AND STORAGE

Store the kit at 2...8 °C. The opened reagents are stable up to the expiry date stated on the label when stored at 2...8 °C.

6. REAGENT PREPARATION

It is very important to bring all reagents and samples to room temperature (20...25 °C) and mix them before starting the test run!

6.1. Microtiterplate

The break-apart snap-off strips are coated with Brucella antigens. Immediately after removal of the strips, the remaining strips should be resealed in the aluminium foil along with the desiccant supplied and stored at 2...8 °C.

6.2. Washing Buffer (20x conc.)

Dilute Washing Buffer 1 + 19; e. g. 10 mL Washing Buffer + 190 mL distilled water. The diluted buffer is stable for 5 days at room temperature (20...25 °C). In case crystals appear in the concentrate, warm up the solution to 37 °C e.g. in a water bath. Mix well before dilution.

6.3. TMB Substrate Solution

The reagent is ready to use and has to be stored at 2...8 °C, away from the light. The solution should be colourless or could have a slight blue tinge. If the substrate turns into blue, it may have become contaminated and should be thrown away.

7. SAMPLE COLLECTION AND PREPARATION

Use serum or milk samples with this assay.

- Serum: Use bovine or porcine serum samples with this assay.
- Pooled serum: Use up to five pooled samples. A maximum of five samples can be pooled.
- Milk: Use bovine milk samples with this assay.

Whole-milk samples can be used after centrifugation for 15 minutes at 2000 x g or 2 minutes at 10.000 x g. The sample should be drawn from below the cream layer.

Avoid repeated freezing and thawing.

Heat inactivation of samples is not recommended.

7.1. Sample Dilution

Before assaying all **serum samples** should be diluted 1+100 with Sample Dilution Buffer. Dispense 10 µL sample and 1 mL Sample Dilution Buffer into tubes to obtain a 1+100 dilution and thoroughly mix with a Vortex.

Pooled serum samples should be diluted with Sample Dilution Buffer.

Dispense e.g. 10 µL of each sample (up to 5 samples) together in 1 mL Sample Dilution Buffer into one tube and thoroughly mix with a Vortex.

Before assaying all **milk samples** should be diluted 1+10 with Sample Dilution Buffer. Dispense 100 µL milk and 1 mL Sample Dilution Buffer into tubes to obtain a 1+10 dilution and thoroughly mix with a Vortex.

Optimal milk dilution might differ from region to region and between animal species. It may range from undiluted to 1+20 and has to be determined by each lab by itself. The described 1+10 dilution has to be considered as a starting point for evaluation.

To the online shop:



8. ASSAY PROCEDURE

Please read the instruction for use carefully **before** performing the assay. Result reliability depends on strict adherence to the instruction for use as described. The following test procedure is only validated for manual procedure. If performing the test on ELISA automatic systems we recommend increasing the washing steps from three up to five and the volume of Washing Buffer from 300 µL to 350 µL to avoid washing effects. Pay attention to chapter 12. Prior to commencing the assay, the distribution and identification plan for all samples and standards/controls (duplicates recommended) should be carefully established on the plate layout supplied in the kit. Select the required number of microtiter strips or wells and insert them into the holder.

Perform all assay steps in the order given and without any delays.

A clean, disposable tip should be used for dispensing each standard/control and sample.

Adjust the incubator to 37 ± 1 °C.

1. Dispense 100 µL standards/controls and diluted samples into their respective wells. Leave well A1 for the Substrate Blank.
2. Cover wells with the foil supplied in the kit.
3. **Incubate for 1 hour ± 5 min at 37 ± 1 °C.**
4. When incubation has been completed, remove the foil, aspirate the content of the wells and wash each well three times with 300 µL of Washing Buffer. Avoid overflows from the reaction wells. The interval between washing and aspiration should be > 5 sec. At the end carefully remove remaining fluid by tapping strips on tissue paper prior to the next step!
Note: Washing is important! Insufficient washing results in poor precision and false results.
5. Dispense 100 µL Conjugate into all wells except for the Substrate Blank well A1.
6. **Incubate for 30 min at room temperature (20...25 °C).** Do not expose to direct sunlight.
7. Repeat step 4.
8. Dispense 100 µL TMB Substrate Solution into all wells.
9. **Incubate for exactly 15 min at room temperature (20...25 °C) in the dark.** A blue colour occurs due to an enzymatic reaction.
10. Dispense 100 µL Stop Solution into all wells in the same order and at the same rate as for the TMB Substrate Solution, thereby a colour change from blue to yellow occurs.
11. Measure the absorbance at 450/620 nm within 30 min after addition of the Stop Solution.

8.1. Measurement

Adjust the ELISA Microtiterplate reader **to zero** using the **Substrate Blank**.

If - due to technical reasons - the ELISA Microtiterplate reader cannot be adjusted to zero using the Substrate Blank, subtract its absorbance value from all other absorbance values measured in order to obtain reliable results!

Measure the absorbance of all wells at **450 nm** and record the absorbance values for each standard/control and sample in the plate layout.

Bichromatic measurement using a reference wavelength of 620 nm is recommended.

Where applicable calculate the **mean absorbance values** of all duplicates.

9. RESULTS

9.1. Run Validation Criteria

In order for an assay to be considered valid, the following criteria must be met:

- **Substrate Blank:** Absorbance value < **0.100**
- **Negative Control:** Absorbance value < **0.200** and < **Cut-off**
- **Cut-off Control:** Absorbance value **0.150 – 1.300**
- **Positive Control:** Absorbance value > **Cut-off**

If these criteria are not met, the test is not valid and must be repeated.

9.2. Calculation of Results

The Cut-off is the mean absorbance value of the Cut-off Control determinations.

Example: Absorbance value Cut-off Control 0.44 + absorbance value Cut-off control 0.42 = 0.86 / 2 = 0.43
Cut-off = 0.43

9.2.1. Results in Units [NTU]

$\frac{\text{Sample (mean) absorbance value} \times 10}{\text{Cut-off}} = [\text{NovaTec Units} = \text{NTU}]$

Example: $\frac{1.591 \times 10}{0.43} = 37 \text{ NTU}$

To the online shop:



9.3. Interpretation of Results

Normal value ranges for this ELISA should be established by each laboratory based on its own sample populations in the geographical areas serviced.

The following values should be considered as a guideline:

Cut-off	10 NTU	-
Positive	> 11 NTU	Antibodies against the pathogen are present. There has been a contact with the antigen (pathogen resp. vaccine).
Equivocal	9 – 11 NTU	Antibodies against the pathogen could not be detected clearly. It is recommended to repeat the test with a fresh sample in 2 to 4 weeks. If the result is equivocal again the sample is judged as negative .
Negative	< 9 NTU	The sample contains no antibodies against the pathogen. A previous contact with the antigen (pathogen resp. vaccine) is unlikely.

Diagnosis of an infectious disease should not be established on the basis of a single test result. A precise diagnosis should take into consideration clinical history, symptomatology as well as serological data.

10. SPECIFIC PERFORMANCE CHARACTERISTICS

The results refer to the groups of samples investigated; these are not guaranteed specifications.

The performance data have been established with selected serum samples. Due to the nature of the Protein A/G conjugate this ELISA should react with other mammalian species also. More detailed information is available on request.

10.1. Precision

Intraassay	n	Mean (E)	CV (%)
#1	24	0.577	4.14
#2	24	1.276	3.34
#3	24	1.200	2.75

Interassay	n	Mean (NTU)	CV (%)
#1	12	23.22	4.97
#2	12	20.13	6.05
#3	12	5.10	8.55

10.2. Diagnostic Specificity

The diagnostic specificity is defined as the probability of the assay of scoring negative in the absence of the specific analyte.

Diagnostic Specificity (bovine and porcine): > 98 % (95 % confidence interval: 88.78 % - 98.62 %)

10.3. Diagnostic Sensitivity

The diagnostic sensitivity is defined as the probability of the assay of scoring positive in the presence of the specific analyte.

Diagnostic Sensitivity (bovine and porcine): > 98 % (95 % confidence interval: 63.37 % - 98.62 %)

10.4. Interferences

Interferences with hemolytic, lipemic or icteric samples are not observed up to a concentration of 10 mg/mL hemoglobin, 5 mg/mL triglycerides and 0.5 mg/mL bilirubin.

10.5. Cross Reactivity

Cross reactions cannot be excluded.

11. LIMITATIONS OF THE PROCEDURE

Bacterial contamination or repeated freeze-thaw cycles of the sample may affect the absorbance values.

To the online shop:



12. PRECAUTIONS AND WARNINGS

- Only for veterinary in-vitro diagnostic use.
- All materials of human or animal origin should be regarded and handled as potentially infectious.
- All components of human origin used for the production of these reagents have been tested for anti-HIV antibodies, anti-HCV antibodies and HBsAg and have been found to be non-reactive.
- Do not interchange reagents or strips of different production lots.
- No reagents of other manufacturers should be used along with reagents of this test kit.
- Do not use reagents after expiry date stated on the label.
- Use only clean pipette tips, dispensers, and lab ware.
- Do not interchange screw caps of reagent vials to avoid cross-contamination.
- Close reagent vials tightly immediately after use to avoid evaporation and microbial contamination.
- After first opening and subsequent storage check conjugate and standard/control vials for microbial contamination prior to further use.
- To avoid cross-contamination and falsely elevated results pipette patient samples and dispense reagents without splashing accurately into the wells.
- The ELISA is only designed for qualified personnel who are familiar with good laboratory practice.

12.1. Safety note for reagents containing hazardous substances

Reagents may contain CMIT/MIT (3:1) or MIT (refer to 4.1)

Therefore, the following hazard and precautionary statements apply.

Warning



H317	May cause an allergic skin reaction.
P261	Avoid breathing spray
P280	Wear protective gloves/protective clothing.
P302+P352	IF ON SKIN: Wash with plenty of soap and water.
P333+P313	If skin irritation or rash occurs: Get medical advice/attention.
P362+P364	Take off contaminated and Wash it before reuse.

Further information can be found in the safety data sheet.

12.2. Disposal Considerations

Residues of chemicals and preparations are generally considered as hazardous waste. The disposal of this kind of waste is regulated through national and regional laws and regulations. Contact your local authorities or waste management companies which will give advice on how to dispose hazardous waste.

13. ORDERING INFORMATION

Prod. No.: BRUVT0050 VetLine Brucella ELISA (96 Determinations)

To the online shop:



BIBLIOGRAPHY / LITERATUR / BIBLIOGRAPHIE / BIBLIOGRAFIA / BIBLIOGRAFÍA / BIBLIOGRAFIA

Wundt W. Krankheiten durch Brucellen. In: Gsell O. und Mohr W. Infektionskrankheiten. Springer Verlag Berlin Heidelberg (1968).

Godfroid J. Brucellosis in wildlife. Rev Sci Tech. 2002 Aug;21(2):277-86.

Salata, R.A., Ravdin, J.I. (1985) Brucella species (Brucellosis). In: Mandell GL, Douglas RG, Bennett JE, EDS. Principles and Practice of Infectious Diseases. New York: John Wiley, 1283-90

Ariza, J., Pellicer, T., Pallarés, R., Foz, A. and Gudiol, F. (1992) Specific antibody profile in human brucellosis. Clin Infect Dis. 1992 Jan;14(1):131-40.

Brouqui P, Raoult D. Endocarditis due to rare and fastidious bacteria. Clin Microbiol Rev. 2001 Jan;14(1):177-207.

Gad El-Rab MO, Kambal AM. Evaluation of a Brucella enzyme immunoassay test (ELISA) in comparison with bacteriological culture and agglutination. J Infect. 1998 Mar;36(2):197-201.

Corbel MJ. Brucellosis: an overview. Emerg Infect Dis. 1997 Apr-Jun;3(2):213-21.





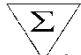
ABBREVIATIONS / ABKÜRZUNGEN / ABRÉVIATIONS / ABBREVIAZIONI / ABREVIACIONES / ABREVIATURAS

CMIT	5-chloro-2-methyl-4-isothiazolin-3-one
MIT	2-methyl-2H-isothiazol-3-one

To the online shop:



SYMBOLS KEY / SYMBOLSCHLÜSSEL / EXPLICATION DES SYMBOLES / LEGENDA / SIMBOLOS / TABELA DE SIMBOLOS

	Manufactured by / Hergestellt von / Fabriqué par / Prodotto da / Fabricado por / Fabricado por
LOT	Lot Number / Chargenbezeichnung / Numéro de lot / Lotto / Número de lote / Número de lote
	Expiration Date / Verfallsdatum / Date de péremption / Scadenza / Fecha de caducidad / Data de Validade
	Storage Temperature / Lagertemperatur / Température de conservation / Temperatura di conservazione / Temperatura de almacenamiento / Temperatura de Armazenamento
REF	Catalogue Number / Katalog Nummer / Référence du catalogue / Numero di codice / Número de Catálogo / Número de Catálogo
	Consult Instructions for Use / Arbeitsanleitung beachten / Consulter la notice d'utilisation / Consultare le istruzioni per l'uso/ Consulte las Instrucciones de Uso / Consultar as Instruções de Utilização
MTP	Microtiterplate / Mikrotiterplatte / Plaque de Microtitrage / Piastre di Microtitolazione / Placa de Microtitulación / Placa de Microtitulação
CONJL	Conjugate / Konjugat / Conjugué / Coniugato / Conjugado / Conjugado
CONTROL -	Negative Control / Negativkontrolle / Contrôle Négatif / Controllo Negativo / Control Negativo / Controle Negativo
CONTROL +	Positive Control / Positivkontrolle / Contrôle Positif / Controllo Positivo / Control Positivo / Controle Positivo
CUT OFF	Cut-off Control / Cut-off Kontrolle / Contrôle Cut-off / Controllo Cut-off / Control Cut-off / Controle Cut-off
DIL	Sample Dilution Buffer / Probenverdünnungspuffer / Tampon de Dilution d'Échantillon / Tampone di Diluizione del Campione / Tampón de Dilución de Muestras / Tampão de Diluição de Amostra
SOLN STOP	Stop Solution / Stopplösung / Solution d'Arrêt / Soluzione Bloccante / Solución de Parada / Solução de Bloqueio
SUB TMB	TMB Substrate Solution / TMB-Substratlösung / Solution de Substrat TMB / Soluzione Substrato TMB / Solución de Sustrato de TMB / Solução Substrato TMB
WASH BUF 20x	Washing Buffer 20x concentrated / Waschpuffer 20x konzentriert / Tampon de Lavage concentré 20 x / Tampone di Lavaggio concentrazione x20 / Tampón de Lavado concentrado x20 / Tampão de Lavagem concentrada 20x
	Contains sufficient for "n" tests / Ausreichend für "n" Tests / Contenu suffisant pour "n" tests / Contenido suficiente per "n" saggi / Contenido suficiente para "n" tests / Conteúdo suficiente para "n" testes

To the online shop:



